



Temperature control from
above and below

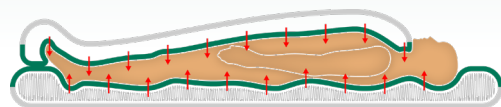

ADULT II (EW II)

Size (L×W) Upper section 70×170cm,
Lower section 115×64cm | Weight approx. 1.540g
Type no. 904

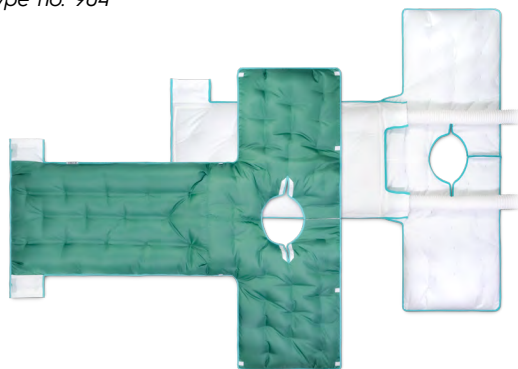

TWINWARM BB – AIR WARMING DEVICE WITH TWO HOSE CONNECTIONS

Temperature setting accurate to the degree between
28 °C and 43 °C | Individual five-stage airflow output
Article no. 12BB-EN-[X]*



TWINWARM

OUR COMMITMENT TO LIFE AND EARTH

Sustainability begins with every decision and does not end with our products. With the MOECK WARMING SYSTEM®, we not only optimize patient care, but also rely on reusable components that combine environmental protection and efficiency.

To this end, we are continuously developing innovative solutions that require fewer resources, promote reuse and enable responsible recycling.

We design our production processes to be as environmentally friendly as possible – for a better future for people and the planet. Because we are convinced that every step counts when it comes to making the world a better place to live.



* X = Please specify the destination country and plug type. An overview of available variants is available on request.
MD002en Rev.02 © 01/2025 - MOECK & MOECK GmbH - All information is without guarantee. We reserve the right to make technical changes to the products as well as printing errors and colour deviations.



MOECK & MOECK

MOECK & MOECK GmbH
Waidmannstraße 12d • 22769 Hamburg • Germany


Tel. +49 (0) 40 4111 4111 • Fax +49 (0) 40 4116 7354
info@moeckundmoeck.de

MOECKWARMINGSYSTEM.COM

MOECK WARMING SYSTEM®

PATIENT CARE ON A HIGH LEVEL

With its comprehensive product solutions, the MOECK WARMING SYSTEM® offers an effective and flexible tool for optimal temperature control of adult and paediatric patients.

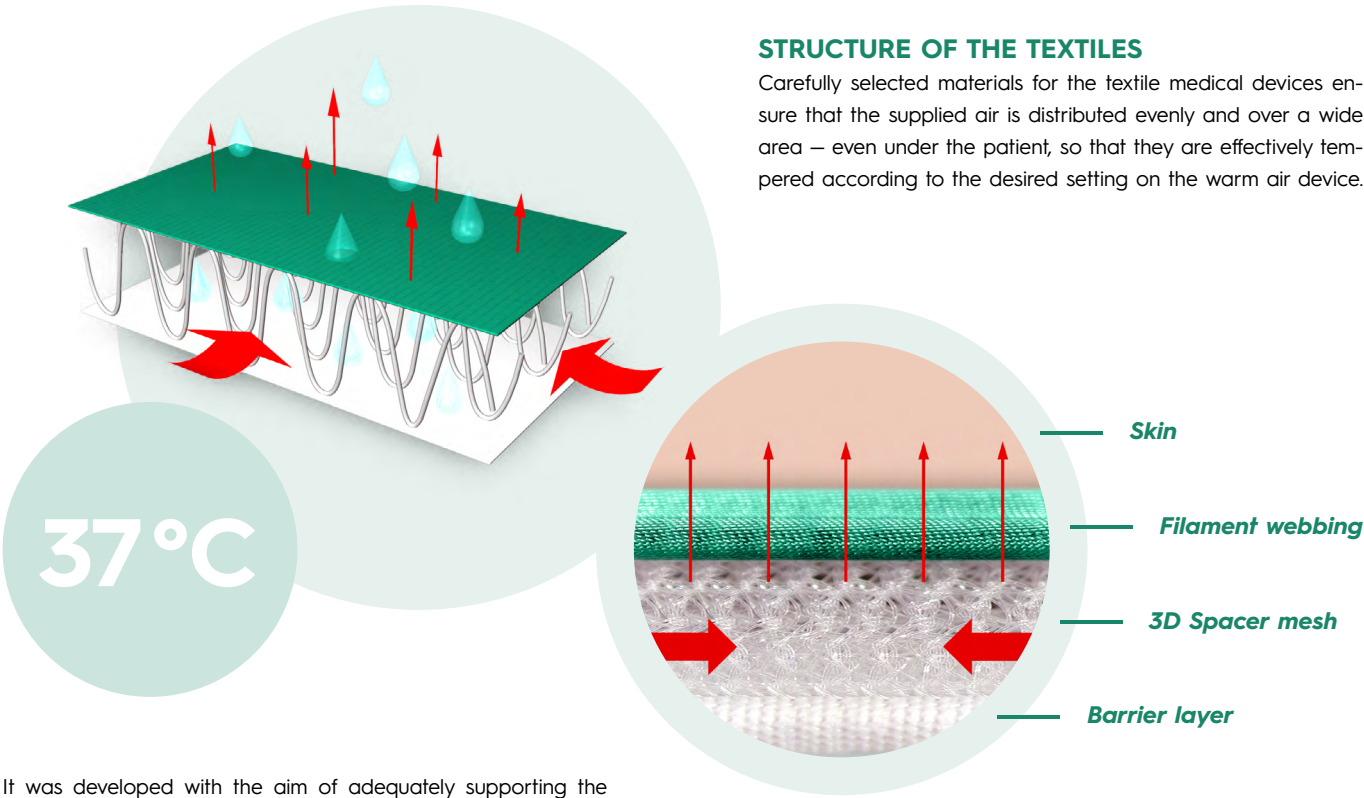

It was developed with the aim of adequately supporting the daily work in medical institutions with a focus on quality-oriented patient care.


SIMPLE BUT EFFECTIVE

The MOECK WARMING SYSTEM® is based on the use of convective air transmission from a air warming device and different models of reusable textile warming blankets as well as warming mats.

During the application, the air tempered by the air warming device is transmitted via one or two flexible hoses into the interior of the respectively connected warming blanket, respectively warming mat.

The integrated layer from the 3D spacer mesh reduces the contact pressure and minimizes fluid accumulation on the surface of the warming mat.

The textile materials used in our warming blankets and warming mats are X-ray-permeable and latex-free. All models are equally suitable for both short and long application periods.

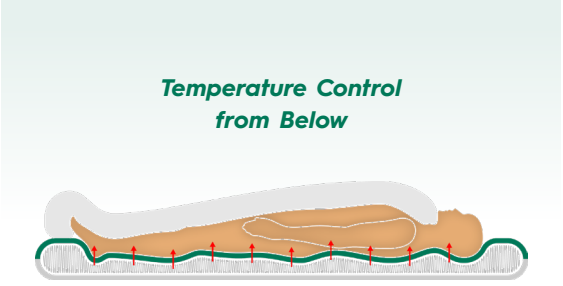


UNIVERSAL II (UNI II)

Size (L×W) 204×64 cm | Weight approx. 1.960g
Type no. 906





UNIVERSAL III (UNI III)

Size (L×W) 115×64 cm | Weight approx. 1.200g
Type no. 910

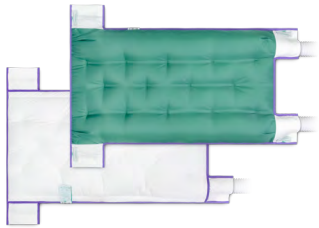


BABY / INFANT (BKK)

Size (L×W) 130×52 cm, Wings 96×49 cm
Weight approx. 890g – Type no. 902

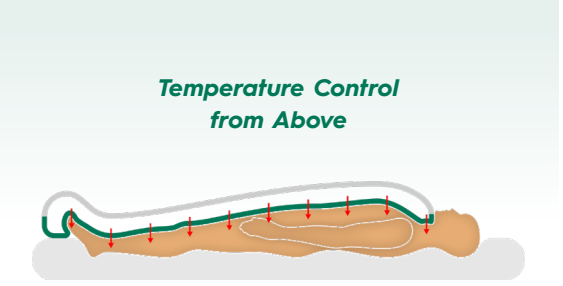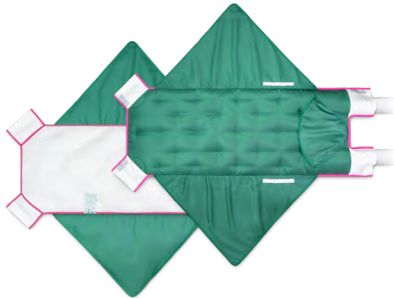


UPPER BODY (OKD)

Size (L×W) 193×70 cm, Head section 57×16 cm
Weight approx. 440g – Type no. 894

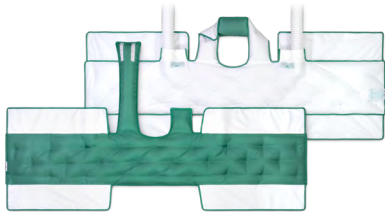

FULL BODY (GKD)

Size (L×W) 196×128 cm | Weight approx. 710g
Type no. 895

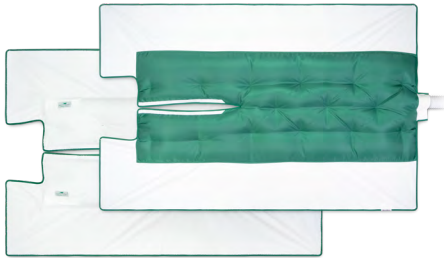

MULTI I

Size (L×W) 180×94 cm, Cut-out 55×36 cm
Weight approx. 520g – Type no. 908

